

# A drug-coated balloon treatment for urethral stricture disease: Two-year results from the ROBUST I study

Rachel A. Mann, MD<sup>1</sup>; Ramón Virasoro, MD<sup>2</sup>; Jessica M. DeLong, MD<sup>2</sup>; Rafael E. Estrella, MD<sup>3</sup>; Merycarla Pichardo, MD<sup>4</sup>; Ramón Rodríguez Lay, MD<sup>5</sup>; Gustavo Espino, MD<sup>6</sup>; Joshua D. Roth, MD<sup>1</sup>; Sean P. Elliott, MD<sup>1</sup>

<sup>1</sup>Department of Urology, University of Minnesota, Minneapolis, MN, United States; <sup>2</sup>Department of Urology, Eastern Virginia Medical School, Norfolk, VA, United States; <sup>3</sup>Clinica Union Medica, Santiago de los Caballeros, Dominican Republic; <sup>4</sup>URUS, Santo Domingo, Dominican Republic; <sup>5</sup>Urology Royal Center Panamá City, Panamá; <sup>6</sup>Centro Especializado San Fernando, Panamá City, Panamá

**Cite as:** Mann RA, Virasoro R, DeLong JM, et al. A drug-coated balloon treatment for urethral stricture disease: Two-year results from the ROBUST I study. *Can Urol Assoc J* 2020 July 27; Epub ahead of print. <http://dx.doi.org/10.5489/auaj.6661>

Published online July 27, 2020

## Abstract

**Introduction:** Mechanical balloon dilation and direct visualization internal urethrotomy (DVIU) are the most widely used treatments for urethral stricture disease in the U.S., but recurrence rates are high, especially after re-treatment. This study investigates the safety and efficacy of the Optilume™ paclitaxel-coated balloon for the treatment of recurrent strictures.

**Methods:** Men with recurrent bulbar strictures  $\leq 2$  cm with 1–4 prior endoscopic treatments were treated with the Optilume™ drug-coated balloon. Patients were evaluated within 14 days, three, six, 12, and 24 months post-treatment. The primary safety endpoint was serious urinary adverse events. The primary efficacy endpoint was  $\geq 50\%$  improvement in International Prostate Symptom Score (IPSS) at 24 months. Secondary outcomes included quality of life, erectile function, flow rate, and post-void residual urine volume.

**Results:** A total of 53 subjects were enrolled and treated; 46 completed the 24-month followup. Forty-three percent of men had undergone  $>1$  previous dilations, with a mean of 1.7 prior dilations. There were no serious adverse events related to treatment at two years. Success was achieved in 32/46 (70%), and baseline IPSS improved from a mean of 25.2 to 6.9 at 24 months ( $p < 0.0001$ ). Quality of life, flow rate, and post-void residual urine volumes improved significantly from baseline. There was no impact on erectile function.

**Conclusions:** Two-year data indicates the Optilume™ paclitaxel-coated balloon is safe for the treatment of recurrent bulbar urethral strictures. Early efficacy results are encouraging and support further followup of these men through five years, as well as further investigation with a randomized trial.

## Introduction

Urethral stricture disease affects approximately 0.6% of males in their lifetime.<sup>1</sup> Direct visualization internal urethrotomy (DVIU) and mechanical dilation remain the most widely used treatments, however, recurrence rates are high when compared to open urethroplasty.<sup>2–4</sup> Moreover, rates of stricture recurrence increase after each endoscopic procedure, making repeat attempts less likely to succeed.<sup>5</sup> Previous studies have investigated the use of anti-proliferative drugs in combination with endoscopic stricture management in an effort to decrease recurrence rates, however, results have been mixed.<sup>6–9</sup> There have also been safety concerns with intralesional injection of mitomycin C (MMC), with reports of serious adverse events.<sup>10</sup>

The Optilume™ drug-coated balloon (DCB; Urotronic, Plymouth, MN) combines urethral dilation with circumferential topical delivery of paclitaxel. Paclitaxel is a microtubule inhibitor with anti-fibrotic and anti-proliferative properties; it is currently used as a coating on vascular stents to prevent restenosis with excellent success.<sup>11–13</sup> Additionally, preliminary animal studies have begun to investigate paclitaxel coating for ureteral stents to prevent stricture after anastomosis.<sup>14</sup> ROBUST I is a single-arm, prospective, multicenter study evaluating outcomes after Optilume™ DCB treatment, with one-year results showing 70% anatomic success with no serious adverse events after 12 months.<sup>15</sup> Herein, we present two-year safety and efficacy outcomes, with efficacy defined as functional success (i.e., symptom score).

## Methods

### Study design and participants

This was a single-arm, prospective, open-label study conducted under a common protocol at four Latin American centers. Eligible patients were men  $\geq 18$  years, with a single

bulbar urethral stricture <12 Fr, and  $\leq 2.0$  cm long on urethrogram. All patients were required to have an International Prostate Symptom Score (IPSS)  $\geq 13$ , a maximum flow rate <10 ml/sec, and to have undergone 1–4 prior endoscopic stricture treatments (none within three months prior to enrollment). Patients who underwent urethroplasty, radical prostatectomy, penile prosthesis or artificial urinary sphincter, or pelvic radiation were excluded, along with patients who underwent intradetrusor onabotulinum toxin A injection within 12 months of study entry. Patients were ineligible if they had urinary stone passage in the previous six months, a diagnosis of lichen sclerosus, chronic kidney disease or serum creatinine >2 mg/dL, neurogenic bladder, bladder or prostate cancer in the previous five years, or active non-genitourinary cancer.

## Procedure

All patients underwent baseline retrograde urethrography prior to study initiation. Prior to DCB treatment, strictures were pre-treated with uncoated balloon dilation and/or DVIU until the urethral lumen diameter increased by 50%. Although the DCB is intended to be used without pre-treatment, pre-treatment was performed in this first-in-man study: 1) to ensure the drug coating was not disrupted when the DCB crossed the stricture; and 2) to prevent double-dosing in patients who had a urethral lumen <20 Fr after the first DCB treatment. The DCB was 3 cm in length; the first several cases used a 24 F balloon until a 30 F became available, after which most cases converted over to that. After pre-treatment, the DCB was inflated to the rated burst pressure and held in place for at least five minutes. Stricture dilation was confirmed fluoroscopically and visually.

## Followup

Patients were seen at five and 14 days, then at three, six, 12, and 24 months with annual followup planned for five years. IPSS, Urethral Stricture Surgery Patient-Reported Outcome Measure (USS-PROM),<sup>16</sup> International Index of Erectile Function (IIEF),<sup>17</sup> maximum flow rate ( $Q_{max}$ ), and post-void residual urine volumes (PVR) were obtained at baseline and at each followup visit. Urethral lumen test (ULT), defined as the ability to pass a flexible cystoscope ( $\geq 15$  Fr) into the bladder or the ability to pass a 14 Fr catheter atraumatically, was performed at six and 12 months post-procedure but not at 24 months; anatomic results by ULT were previously reported. Herein, we focus on 24-month functional outcomes.

## Primary safety endpoint

The primary safety endpoint was the rate of treatment-related urinary severe adverse events (SAEs), defined as urethral

fistula formation, de novo urinary retention >14 days post-treatment, de novo stress incontinence (>1 pad/day), or urethral rupture.

## Efficacy endpoint

The primary efficacy endpoint was defined as a 50% improvement in IPSS compared to baseline in the absence of retreatment. Because symptom scores can fluctuate over time, we report the IPSS score at multiple time points but define success vs. failure on the date of last followup. Patients were considered to have failed treatment if IPSS improvement was not  $\geq 50\%$  at the time of last followup or if the subject required retreatment. If a patient exited the study without  $\geq 50\%$  IPSS improvement but ULT at the time of study exit demonstrated a patent urethra, the patient was not considered a failure, but instead censored, as symptoms were attributable to a different cause, such as prostatic enlargement.

## Secondary endpoints

Secondary endpoints included: 1) improvement in lower urinary tract symptoms (LUTS) based on the USS-PROM, a validated instrument developed to assess urinary symptoms following stricture surgery; 2) sexual function, using the “overall satisfaction” question of the IIEF; and 3)  $Q_{max}$  and PVR.

## Data analysis

Baseline characteristics and the primary and secondary endpoints were tabulated using descriptive analysis. The number and percentage of subjects experiencing at least one device-related SAE were presented, and the number and percentage of IPSS failures, ULT failures, patient requiring retreatment, as well as the change in IPSS from baseline were reported. Additional analysis was performed, comparing one- and two-year success rates based on stricture length, balloon diameter, and number of prior endoscopic treatments. A Chi-squared test was used, with  $p < 0.05$  indicating a statistically significant difference. For secondary endpoints, the mean and standard deviation were presented at various timepoints and compared to baseline values. Patients were censored from the secondary analysis if they underwent retreatment after six months.

## Results

### Patients

Between November 29, 2016 and September 9, 2017, 53 patients were enrolled and treated with the DCB. The aver-

age was 50.7 years (range 22–81), the majority (83%) were Hispanic/Latino, and the mean number of endoscopic treatments prior to enrollment was 1.7. Some patients performed self-dilation independently, however, only physician-administered procedures were recorded. All strictures were located in the bulbar urethra, with an average length of 0.9 cm and diameter of 2.3 mm. Prior to DCB treatment, pre-dilation was accomplished with an uncoated balloon in 59%, DVIU in 15%, or a combination of the two in 26% (Table 1).

### Primary safety and efficacy endpoints

There were no treatment-related urinary SAEs at two years post-procedure. There were six non-urinary SAEs but none were directly attributable to the device or procedure. In total, there were 71 adverse events, most commonly urinary tract infection (17%), fever (8%), dysuria (7%), acute urinary retention (6%), and headache (6%). The majority were classified as mild (61%) or moderate (31%) according to the Common Terminology Criteria for Adverse Events, and 23% were categorized as “possibly,” “probably,” or “definitely” related to the procedure or device (Fig. 1).

Functional treatment success (defined as IPSS improve-

ment  $\geq 50\%$  without retreatment) was achieved in 32/46 (70%) of men at 24 months (Fig. 2). The 14 failures included seven men with IPSS that did not improve  $\geq 50\%$  at last followup and seven men who underwent retreatment (Fig. 3). Five of 53 men were censored due to being lost to followup before two years but with IPSS improvement  $\geq 50\%$  without retreatment and with an open urethra on ULT at the time of exit. One man did not have 50% improvement in IPSS but did have a successful ULT at his six-month visit; he was censored at six months due to a clinical determination that symptoms were due to prostatic enlargement, in the absence of stricture recurrence. One man exited the study within one month of treatment without any data. Balloon size was a statistically significant predictor of treatment success and notably, 11 of the 14 failures at two years occurred with the 24 F balloon ( $p=0.010$ ) (Table 2). Stricture length and number of prior endoscopic treatments did not predict one- or two-year success (Table 2). The mean IPSS among all patients improved from a baseline mean of 25.2 (range 15.0–34.0) to 6.9 (range 0.0–34.0) at 24 months ( $p<0.0001$ ). Additionally, mean IPSS quality of life score also significantly improved from a baseline of 4.9 (range 2.0–6.0) to 0.9 (range 0.0–6.0) (Table 3).

### Secondary endpoints


There were seven patients with anatomic failure based on ULT (six at six months, and one at 12 months). Of the seven patients with anatomic recurrence, three were functional failures and three were functional successes at 24 months; one was censored at six months with IPSS=3. The three functional failures included two who underwent retreatment and one who exited the study without 50% IPSS improvement. Three men had anatomic recurrence but  $\geq 50\%$  IPSS improvement at 24 months (so considered as functional successes). There were 38 patients with a successful ULT; 34 of these completed the 24-month followup. There were eight patients who did not complete the ULT at six or 12 months. Among the ULT successes, mean IPSS improved 19 points compared to baseline (25.8 $\pm$ 4.2 at baseline to 6.8 $\pm$ 8.0 at 24 months). There were five IPSS failures at 24 months among patients who had passed the ULT at 12 months.

Mean USS-PROM scores improved from 15.9 ( $\pm 4.69$ ) at baseline to 3.6 ( $\pm 5.79$ ) at 24 months ( $p<0.05$ ; Table 3). Mean IIEF satisfaction scores did not significantly change: 6.5 ( $\pm 2.6$ ) at baseline and 7.6 ( $\pm 2.48$ ) at 24 months ( $p>0.05$ ) (Table 3). Mean  $Q_{max}$  improved from 5.0 ml/sec (baseline) to 22.2, 20.1, and 17.5 ml/sec at three, 12, and 24 months' followup, respectively ( $p<0.05$ ) (Table 3). Additionally, the mean PVR decreased from 141.0 ml (baseline) to 36.5, 24.6, and 45.5 ml after three, 12, and 24 months ( $p<0.05$ ) (Table 3).

**Table 1. Patient demographics and urologic medical history**

Variables	N=53
Patient demographics	
Age (years)	
Mean $\pm$ SD	50.7 $\pm$ 15.47
Range	22.0–81.0
Median	50.0
Male gender, n (%)	53 (100.0%)
Race of subjects, n (%)	
Black or African origin	8 (15.1%)
Hispanic or Latino	44 (83.0%)
Other	1 (1.9%)
Suprapubic catheter at baseline	7 (13.2%)
Stricture etiology, n (%)	
Iatrogenic	24 (45.3%)
Idiopathic	2 (3.8%)
Traumatic	27 (50.9%)
Stricture measurements, mean $\pm$ SD (n)	
Stricture length (mm)	9.00 $\pm$ 5.20 (53)
Urethral diameter at stricture (mm)	2.47 $\pm$ 1.97 (53)
Urethral diameter at area healthy tissue (mm)	10.20 $\pm$ 3.62 (53)
Pre-treatment	
Uncoated balloon	31/53 (59%)
DVIU	8/53 (15%)
Uncoated balloon + DVIU	14/53 (26%)
Number of previous endoscopic treatments, n (%)	
1	30 (57%)
2	13 (25%)
3	9 (17%)
4	1 (2%)

DVIU: direct visualization internal urethrotomy; SD: standard deviation.

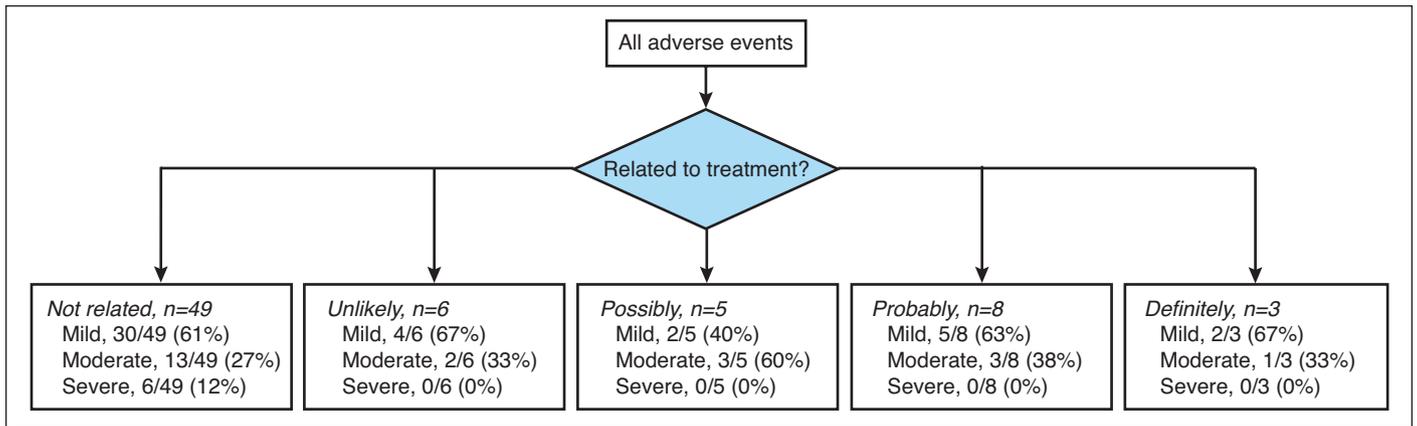

**Fig. 1.** Adverse events, relation to treatment, and severity. Severity was determined using the Common Terminology Criteria or Adverse Events (CTCAW, v4.0) Severity Scale.


### Discussion

The ROBUST I trial is a multicenter, single-arm clinical trial investigating the safety and efficacy of the Optilume™ DCB among male patients with a single, recurrent, bulbar urethral stricture. At two years, there were no treatment-related SAEs and the 24-month success rate was 70%, defined as IPSS improvement  $\geq 50\%$  in the absence of retreatment.


Treatment safety is of paramount concern, especially in light of reports of SAEs following intraurethral MMC injection, including urethral fistula formation, osteitis pubis, and tissue necrosis.<sup>10</sup> In our study, the majority of adverse events were mild to moderate in severity, and the most common treatment-related adverse event was urinary tract infection. The Optilume™ DCB may have a lower risk of SAE compared to intralesional injections because the drug is applied topically; the mechanical balloon dilation ruptures the urethral mucosa allowing direct drug delivery to the exposed submucosa through fissures in the stricture. This controlled drug delivery system limits dose variability between users.

Additionally, Good Laboratory Practice Animal study data acquired prior to study initiation revealed that urethral tissue concentration of paclitaxel dropped 73% after seven days.<sup>18</sup> Serum drug levels in the 53 men in this study were also very low at all times points.

At two years, stricture recurrence rate based on IPSS or retreatment remained low (30%). Mean IPSS improved significantly from 25.2 prior to DCB treatment to 6.9 at 24 months. Patients who underwent dilation with the 30 F balloon were significantly more likely to have continued treatment success at two years. Based on these results, we are currently recommending only the 30 F balloon for bul-

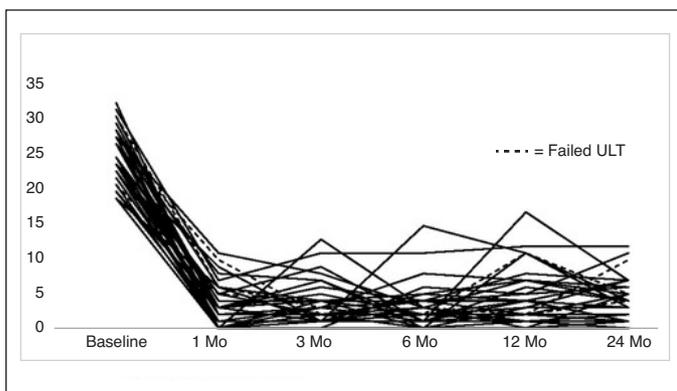

**Fig. 2.** IPSS over time among treatment successes, defined as IPSS improvement  $\geq 50\%$  after 24 months in the absence of re-treatment. Three patients failed the ULT at 6 or 12 months but still met criteria for functional success at 24 months. ---: failed ULT. IPSS: International Prostate Symptom Score; ULT: urethral lumen test.

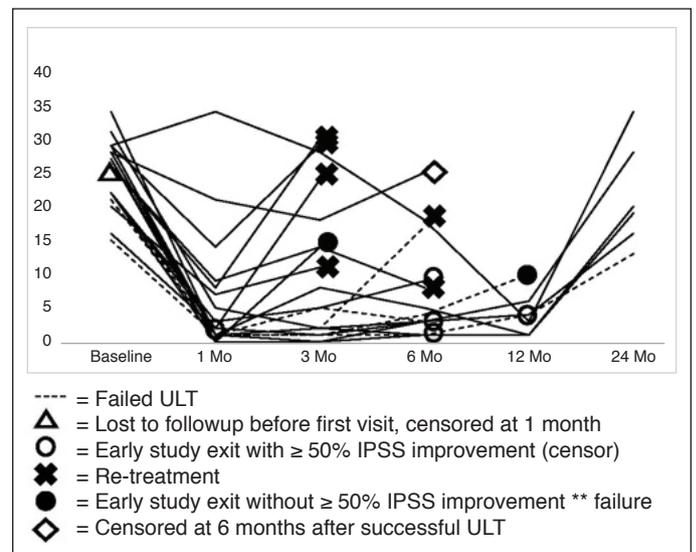

**Fig. 3.** IPSS over time among censored patients and treatment failures, defined as failure to improve IPSS improvement  $\geq 50\%$  at the time of study exit or re-treatment. There were four patients who failed the ULT in this cohort; one patient underwent re-treatment, one patient was censored at six months with an IPSS of 1, and two patients exited the study without  $\geq 50\%$  IPSS improvement. All patients were considered failures at 24 months based on failure to improve IPSS  $\geq 50\%$ . IPSS: International Prostate Symptom Score; ULT: urethral lumen test.

**Table 2. Device and patient characteristics related to functional success**

	2-year success, n/N (%)	p
Drug-coated balloon size		0.010*
8 mm (n=25)	12/23 (52%)	
10 mm (n=28)	20/23 (87%)	
Stricture length		0.386
≤0.5 cm (n=16)	6/11 (55%)	
0.6–1.0 cm (n=22)	14/21 (67%)	
1.1–1.5 cm (n=9)	7/8 (88%)	
1.6–2.0 cm (n=6)	5/6 (83%)	
Number of prior endoscopic treatments		0.090
1 (n=29)	20/25 (80%)	
2 (n=14)	6/13 (46%)	
≥3 (n=10)	6/8 (75%)	

p was obtained using a Chi-squared analysis. \*Denotes statistical significance (p≤0.05).

bar strictures in the ROBUST III randomized trial; we reserve the 24 F balloon for the penile urethra. Secondary outcome measures demonstrated a durable improvement in USS-PROM scores and PVR from baseline.  $Q_{max}$  improved from baseline (5 mL/sec) at all time points but generally decreased from the first post-procedure measure at three months (22 mL/sec) to the last measure at 24 months (18 mL/sec), despite IPSS scores remaining low at all postoperative time points. We will follow this trend in  $Q_{max}$  through five years. Erectile function appeared to be unaffected by treatment (Table 3).

Our study is designed to conduct annual followup for five years, and we now present two-year preliminary data. One-year efficacy endpoint from the ROBUST I study was defined as anatomic success based on the ULT, regardless of symptoms or flow rate.<sup>15</sup> Patient symptoms are an important measure of treatment success, as anatomic recurrence does not always correlate with worsening LUTS; because of this and in order to avoid invasive procedures, our protocol does not call for performing a ULT after one year and defines long-term success based on objective (retreatment) and subjective (symptomatic) components, similar to much of the urethroplasty literature. In a large, multi-institutional study evaluating outcomes after bulbar urethroplasty, up to 35% of patients with anatomic recurrence on cystoscopy were asymptomatic.<sup>19</sup> Additionally, recent data from the TURNS

study group showed that although cystoscopic stricture recurrence after urethroplasty was associated with secondary interventions, it was a poor predictor of patient symptoms.<sup>20</sup> In our study, four of the seven patients with ULT failure were functional successes at the time of study exit. Recent literature suggests an important role for patient-reported outcome measures; globally, there has been a paradigm shift to place more emphasis on patient symptoms and satisfaction when defining treatment success.<sup>21</sup>

In this early-phase study, one of the notable limitations was the lack of a control arm. We also excluded patients with a history of radiation therapy, penile urethral strictures, bladder neck contractures, or patients with a history of lichen sclerosus; it is unclear how the DCB would have performed in these populations. Urethral lumen test was not performed at the two-year mark, and in some patients, worsening of symptoms may have been due to other causes (i.e., benign prostatic hyperplasia) rather than stricture recurrence. This study is intended to continue for five years, and the 24-month followup does not capture delayed stricture recurrence — although data suggest that most recurrences occur within one year of treatment.<sup>22</sup> Additionally, patients in our cohort performed self-dilation independently prior to study entry and it is uncertain how this may have affected recurrence rates. Lastly, the DCB is intended to be used without pre-dilation, and future studies will be necessary to determine if the efficacy of the DCB is affected when used alone.

## Conclusions

The Optilume™ DCB is safe and symptomatic evidence of stricture recurrence is low at two years post-procedure. Data from the ROBUST I study will be reported through five years and a randomized clinical trial comparing the DCB to DVIU or plain balloon dilation is ongoing.

**Competing interests:** Dr. Virasoro, Dr. DeLong, and Dr. Elliott serve as consultants for Urotronic. The remaining authors report no competing personal or financial interests related to this work.

**Table 3. Primary and secondary outcome measures over time**

Category	Baseline	3 months	6 months	12 months	24 months
IPSS, mean ± SD (n)	25.2±4.46 (53)	6.1±7.63 (51)	4.6±5.18 (45)	4.5±3.90 (40)	6.9±7.66* (38)
IPSS QOL, mean ± SD (n)	4.9±0.86 (53)	0.8±1.32 (51)	0.7±0.93 (45)	0.7±0.93 (40)	0.9±1.47* (38)
USS-PROM score, mean ±SD (n)	15.9±4.69 (53)	3.2±5.53 (51)	1.9±2.91 (45)	1.4±1.78 (40)	3.6±5.79* (38)
IIEF: Overall satisfaction					
Mean ±SD (n)	6.5 ± 2.62 (53)	7.9 ± 2.53 (51)	7.6 ±2.80 (45)	8.1±2.51 (40)	7.6±2.48 (38)
Qmax (mL/sec), mean ± SD (n)	5.0±2.56 (46)	22.2±12.54 (51)	19.8±10.78 (45)	20.1±10.0 (39)	17.5±10.39* (38)
PVR (mL), mean ± SD (n)	141.4±105.05 (43)	36.5±37.66 (51)	30.0±42.84 (45)	24.6±32.08 (39)	45.5± 49.50* (38)

Patients who underwent retreatment after 6 months were censored from the analysis. \*Statistically significant change from baseline value (p≤0.05). IIEF: International Index of Erectile Function; IPSS: International Prostate Symptom Score; PVR: post-void residual; Qmax: maximum flow rate; SD: standard deviation; ULT: urethral lumen test; USS-PROM: Urethral Stricture Surgery Patient-Reported Outcome Measure.

**Funding:** This study was supported by Urotronic.

This paper has been peer-reviewed

## References

1. Santucci RA, Joyce GF, Wise M. Male urethral stricture disease. *J Urol* 2007;177:1667-74. <https://doi.org/10.1016/j.juro.2007.01.041>
2. Liu JS, Hofer MD, Oberlin DT, et al. Practice patterns in the treatment of urethral stricture among American urologists: A paradigm change? *Urology* 2015;86:830-4. <https://doi.org/10.1016/j.urology.2015.07.020>
3. Meeks JJ, Erickson BA, Granieri MA, et al. Stricture recurrence after urethroplasty: a systemic review. *J Urol* 2009;182:1266-70. <https://doi.org/10.1016/j.juro.2009.06.027>
4. Shaw NM, Venkatesan K. Endoscopic management of urethral stricture: Review and practice algorithm for management of male urethral stricture disease. *Curr Urol Rep* 2018;19:19. <https://doi.org/10.1007/s11934-018-0771-6>
5. Steenkamp JW, Heyns CF, de Kock ML. Internal urethrotomy versus dilation as treatment for male urethral strictures: A prospective, randomized comparison. *J Urol* 1997;157:98-101. [https://doi.org/10.1016/S0022-5347\(01\)65296-0](https://doi.org/10.1016/S0022-5347(01)65296-0)
6. Farrell MR, Lawrenz CW, Levine LA. Internal urethrotomy with intralesional mitomycin C: An effective option for endoscopic management of recurrent bulbar and bulbomembranous urethral strictures. *Urology* 2017;110:S223. <https://doi.org/10.1016/j.urology.2017.07.017>
7. Mazdak H, Izadpanahi MH, Ghalamka A, et al. Internal urethrotomy and intravaginal submucosal injection of triamcinolone in short bulbar urethral strictures. *Int Urol Nephrol* 2010;41:565-8. <https://doi.org/10.1007/s11255-009-9663-5>
8. Tavakkoli Tabassi K, Yarmohamadi A, et al. Triamcinolone injection following internal urethrotomy for treatment of urethral stricture. *Urol J* 2011;8:132-6.
9. Modh R, Cai PY, Sheffield A, et al. Outcomes of direct vision internal urethrotomy for bulbar urethral strictures: Technique modification with high-dose triamcinolone injection. *Adv Urol* 2015;2015:281969. <https://doi.org/10.1155/2015/281969>
10. Intralesional injection of mitomycin C at transurethral incision of bladder neck contracture may offer limited benefit: TURNS Study Group. *J Urol* 2015;193:587-92. <https://doi.org/10.1016/j.juro.2014.08.104>
11. Käsmann L, Manig L, Janssen S, et al. Chemoradiation including paclitaxel for locally recurrent muscle-invasive bladder cancer in elderly patients. *In Vivo* 2017;31:239-42. <https://doi.org/10.21873/invivo.11051>
12. Herten M, Torsello GB, Stahlhoff S. Critical appraisal of paclitaxel balloon angioplasty for femoral-popliteal arterial disease. *Vasc Health Risk Manag* 2016;12:341-56. <https://doi.org/10.2147/VHRM.S81122>
13. Habib A, Finn AV. Antiproliferative drugs for restenosis prevention. *Interv Cardiol Clin* 2016;5:321-9. <https://doi.org/10.1016/j.iccl.2016.02.002>
14. Kram W, Rebl H, Wyrwa R, et al. Paclitaxel-coated stents to prevent hyperplastic proliferation of ureteral tissue: From in vitro to in vivo. *Urolithiasis* 2020;48:47-56. <https://doi.org/10.1007/s00240-018-1081-7>
15. Virasoro R, DeLong JM, Mann RA, et al. A drug-coated balloon treatment for urethral stricture disease: Interim results from the ROBUST I study. *Can Urol Assoc J* 2020;14:187-91. <https://doi.org/10.5489/cuaj.6323>
16. Jackson MJ, Sciberras J, Magera A, et al. Defining a patient-reported outcome measure for urethral stricture surgery. *Eur Urol* 2011;60:60-8. <https://doi.org/10.1016/j.eururo.2011.03.003>
17. Rosen RC, Riley A, Wagner G, et al. The international index of erectile function (IIEF): A multidimensional scale for assessment of erectile dysfunction. *Urology* 1997;49:822-30. [https://doi.org/10.1016/S0090-4295\(97\)00238-0](https://doi.org/10.1016/S0090-4295(97)00238-0)
18. Barnett P, Virasoro R, Rousselle D. Tissue and pharmacokinetic response to treatment with a drug-coated balloon in a canine urethra model. Poster presented at: World Congress of Endourology, Sep 20, 2018; Paris, France.
19. Erickson BA, Elliott SP, Voelzke BB, et al. Multi-institutional 1-year bulbar urethroplasty outcomes using a standardized prospective cystoscopic followup protocol. *Urology* 2014;84:213-6. <https://doi.org/10.1016/j.urology.2014.01.054>
20. Baradaran N, Fergus KB, Moses RA, et al. Clinical significance of cystoscopic urethral stricture recurrence after anterior urethroplasty: A multi-institution analysis from Trauma and Urologic Reconstructive Network of Surgeons (TURNS). *World J Urol* 2019;37:1263-2768. <https://doi.org/10.1007/s00345-019-02653-6>
21. Baradaran N, Hampson LA, Edwards TC, et al. Patient-reported outcome measures in urethral reconstruction. *Curr Urol Rep* 2018;19(7):48. <https://doi.org/10.1007/s11934-018-0797-9>
22. Heyns CF, Steenkamp JW, Kock ML, et al. Treatment of male urethral strictures: Is repeated dilation or internal urethrotomy useful? *J Urol* 1998;160:356-8. [https://doi.org/10.1016/S0022-5347\(01\)62894-5](https://doi.org/10.1016/S0022-5347(01)62894-5)

**Correspondence:** Dr. Rachel A. Mann, Department of Urology, University of Minnesota, Minneapolis, MN, United States; [mann0469@umn.edu](mailto:mann0469@umn.edu)