


## injeTAK<sup>®</sup>

Cystoscopy needle for precise placement.

### UROLOGY & UROGYNECOLOGY

For cystoscopy injections into the bladder wall for Overactive Bladder (OAB) and Neurogenic Detrusor Overactivity (NDO)

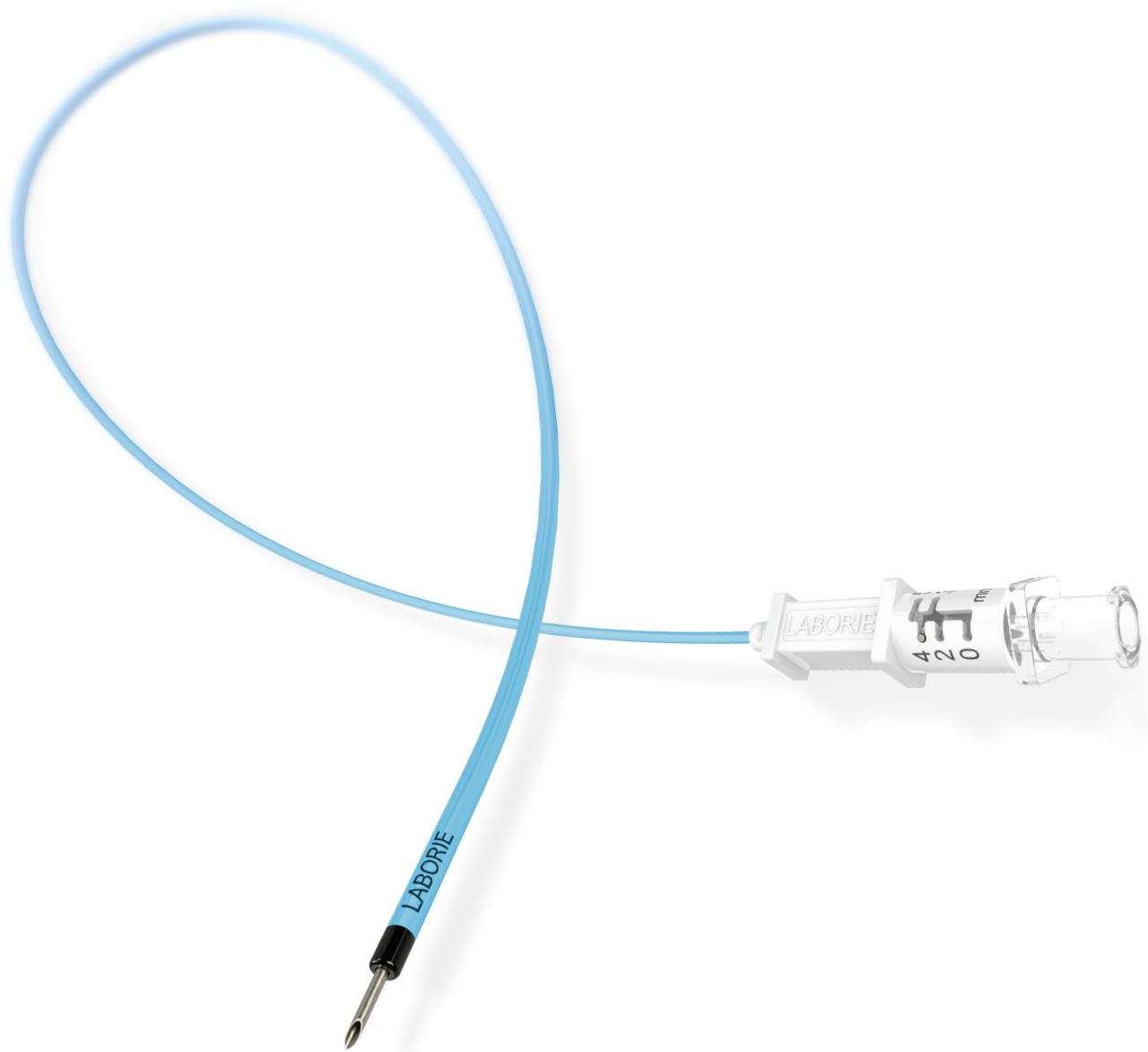

# injeTAK

## Patented depth dialing for accuracy, safety and speed

### Unmatched patient treatment and safety

- For cystoscopy injections into the bladder wall for Overactive Bladder (OAB) and Neurogenic Detrusor Overactivity (NDO)
- Accurate, safe dispensing with decreased procedure time

### Optimized adjustability



- Retractable needle protects needle tip and cystoscope
- Does not require introducer sheath
- Needle tip can be adjusted to 0, 2, 3, 4 and 5 mm to accommodate bladder wall thickness
- Settings provide optimal tip length for speed, safety and consistency (2 mm per RX manufacture's PI)
- Ergonomic design allows for easy tip adjustments

### Needle characteristics

- 4.8 Fr allows for use with smaller working channels and most pediatric scopes
- Black indicator tip provides visual proximity to tissue for accurate placement, safe penetration and gauged puncture depth
- Precise angle and 3-bevel tip engineered to minimize puncture force
- Unique cannula design provides easy injection and superior flexibility while maintaining structural integrity

### Available size options

- 35 cm – DIS199
- 70 cm – DIS201

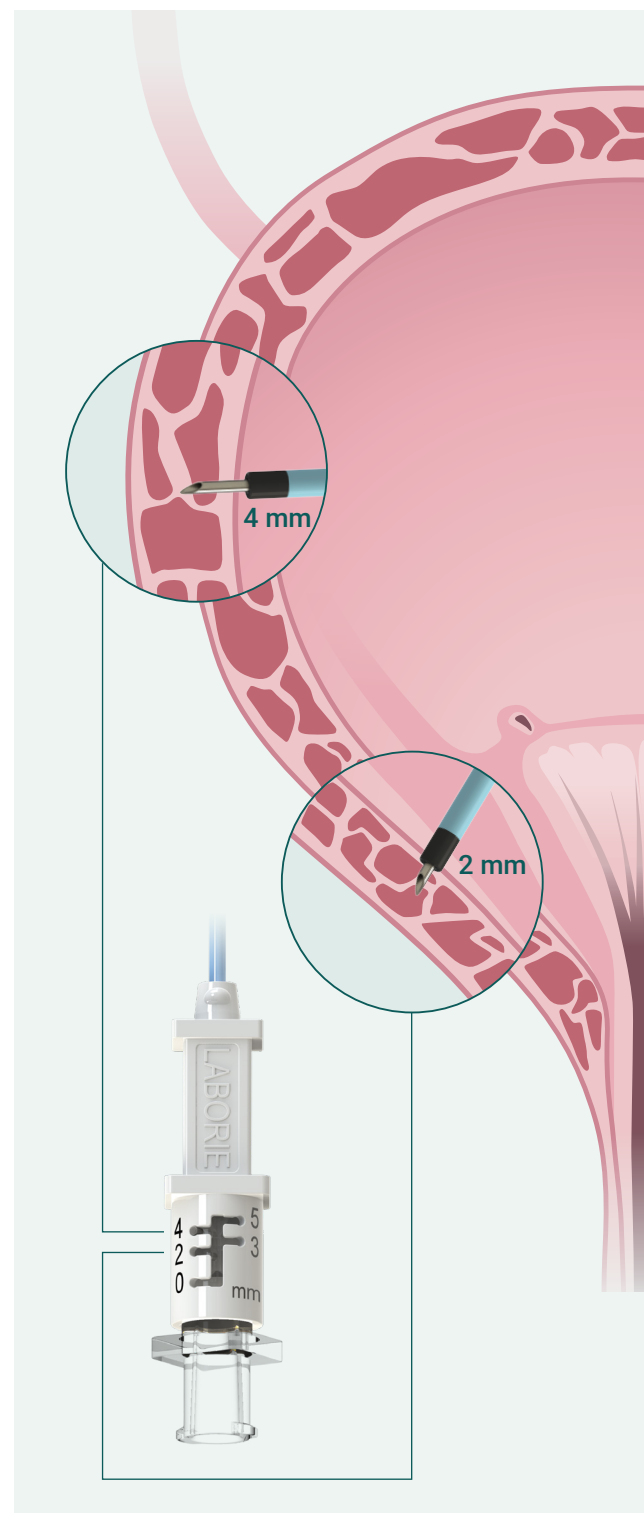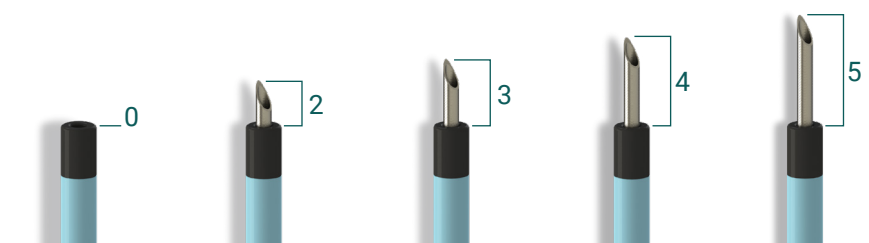

All trademarks are the property of their respective owners. © 2022 Laborie. All rights reserved.

## Laborie<sup>1</sup>

USA

Pease International Tradeport  
180 International Drive,  
Portsmouth, NH 03801, USA

T +1 802 857 1300  
F +1 802 878 1122  
W [laborie.com](http://laborie.com)  
E [usmarketing@laborie.com](mailto:usmarketing@laborie.com)



Learn more

112776-0 10/22