


## PyloPlus

PyloPlus Rapid Urease Test  
Helicobacter Pylori Test

GASTROENTEROLOGY  
Rapid Results in 1 Hour  
No Refrigeration and No Reagents  
3 Year Shelf Life

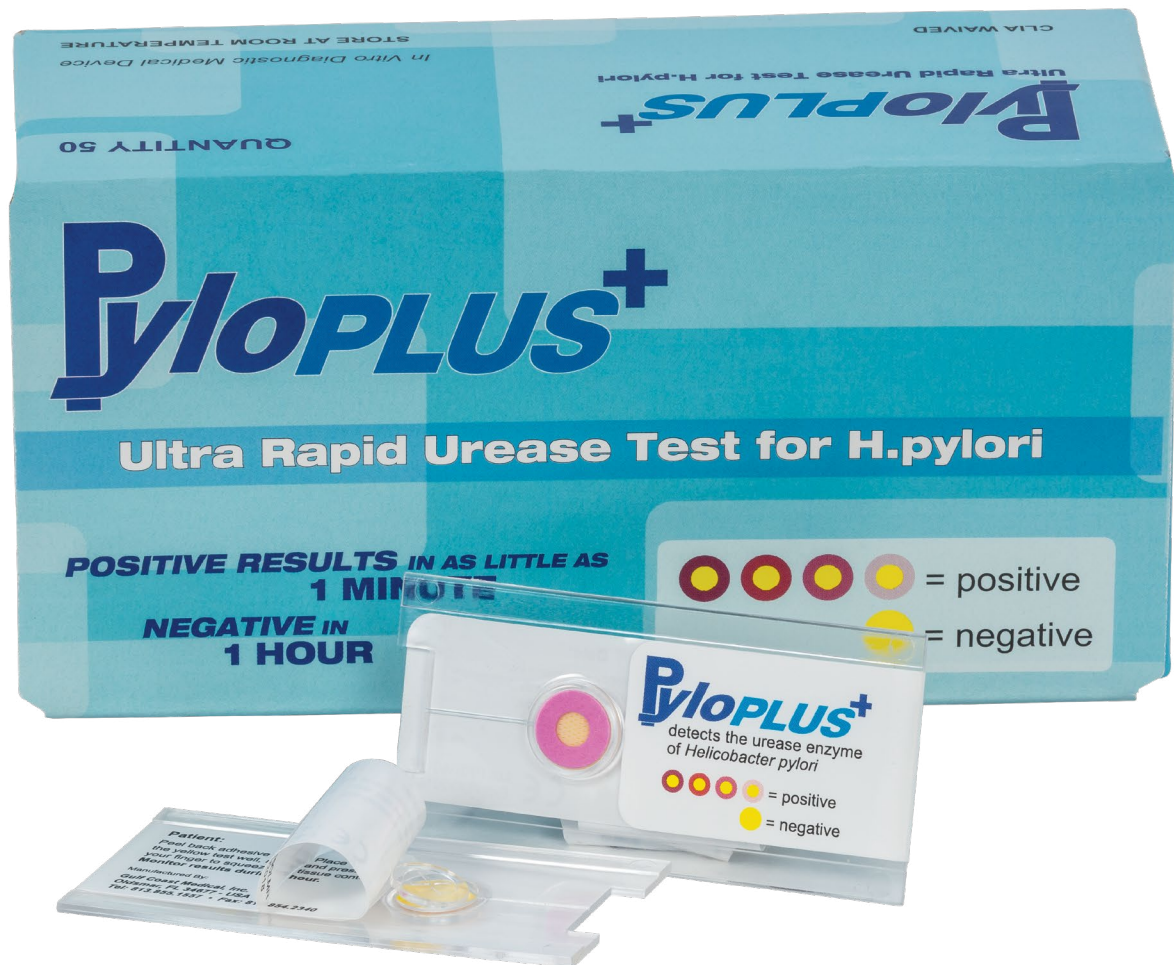

**PyloPLUS<sup>+</sup>**

# PyloPlus: The Newest Generation


## Helicobacter Pylori Test

While there are several methods of detecting H.pylori, biopsy testing using rapid urease is considered the "Gold Standard" by professionals worldwide.

PyloPlus Rapid Urease Test (RUT) is the most effective and easy to use Rapid Urease Test on the market today. PyloPlus RUT offers significant improvements that make detecting H.pylori infections in gastric mucosal biopsies faster, easier, less expensive and more accurate than ever before.




Feature	Benefit
No Refrigeration Required	Can be stored at room temperature. No need to keep temperature logs
Rapid Results	Positive results in as little as 1 minute, negative results in 1 hour
Accurate	99% sensitivity and 97% specificity as shown in major US studies comparing PyloPlus to histology and other rapid urease tests
Dry Test	No incubation or reagents required
Cost Effective	Low per test cost
Simple and Convenient	CLIA waived
3 Year Shelf Life	Never worry about product dating

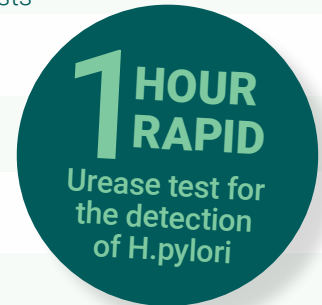


### 5 Simple Steps


1. Peel back label

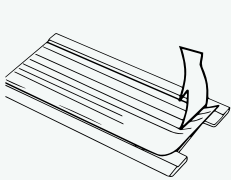

2. Place one or more biopsies on test

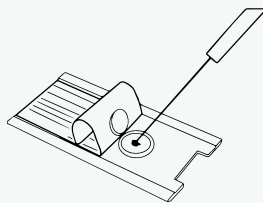

3. Re-seal test

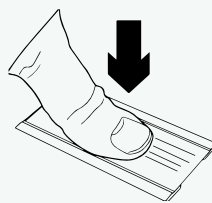

4. Wait up to one hour



5. Read results

Red = Positive



Yellow = Negative



Product No.	Description	Packaging
PYLO-50	PyloPlus Rapid Urease Test	50 each

Laborie is an authorized distributor of PyloPlus for ARJ Medical, Inc © 2022 Laborie Medical Technologies Corp. All Rights Reserved.

**Laborie**

USA

T +1 800 522 6743

E [gis-gisupplyinfo@laborie.com](mailto:gis-gisupplyinfo@laborie.com)  
W [laborie.com](http://laborie.com)



Learn more

109352-1 08/22