


## injeTAK<sup>™</sup>

Cystoscopy needle for precise placement

UROLOGY & UROGYNECOLOGY

Adjustable tip-length cystoscopy  
needles for precision injection into  
the bladder wall

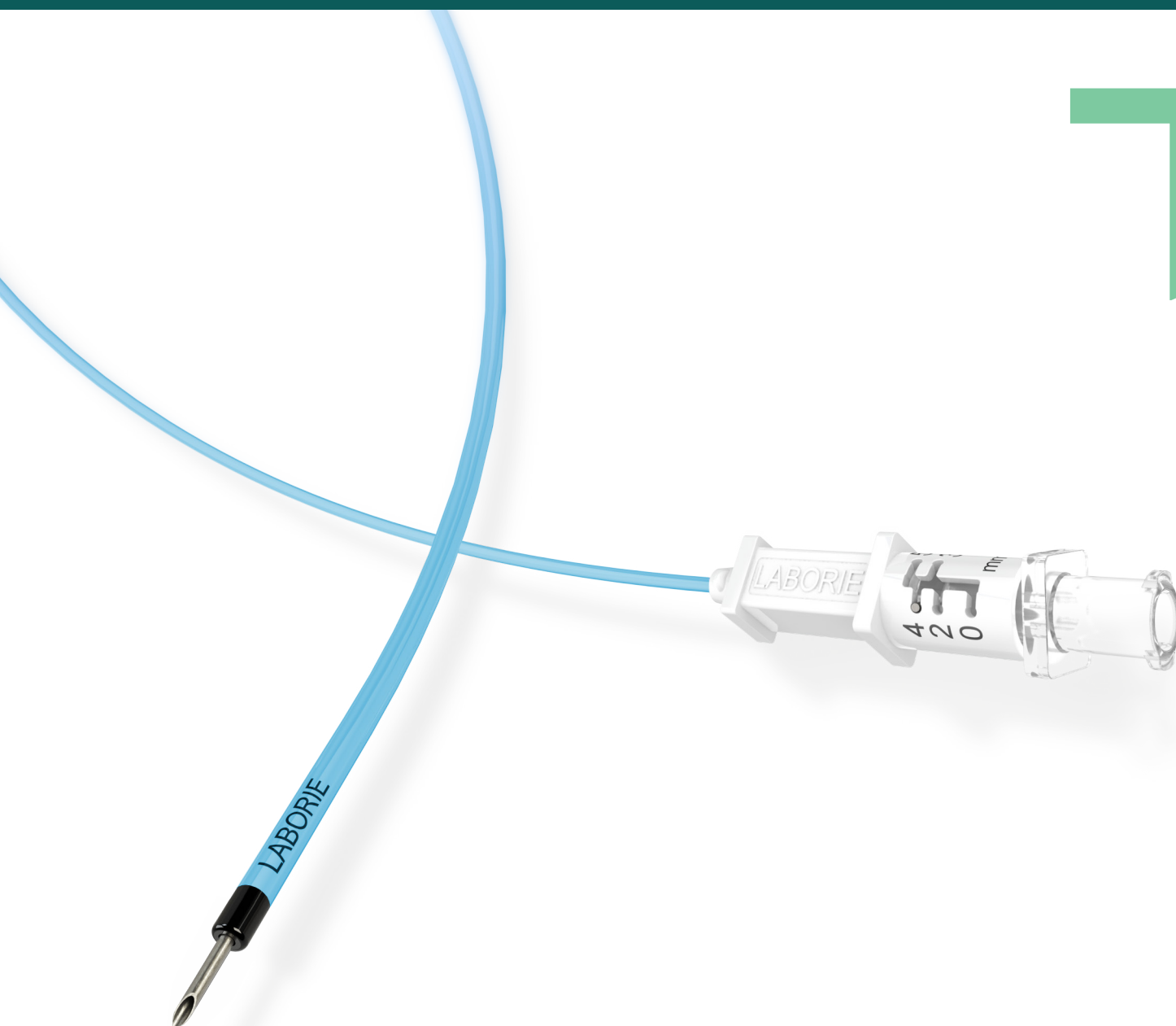

# injeTAK™, patented depth dialing for accuracy and safety

## Unmatched patient treatment and safety

- For cystoscopy injections into the bladder wall for Overactive Bladder (OAB) and Neurogenic Detrusor Overactivity (NDO)
- Accurate and safe dispensing

## Optimized adjustability


- Retractable needle protects needle tip and cystoscope
- Does not require introducer sheath
- Needle tip can be adjusted to 0, 2, 3, 4 and 5 mm to accommodate bladder wall thickness
- Settings provide optimal tip length for safety and consistency (2 mm per RX manufacture's PI)
- Ergonomic design allows for easy tip adjustments

## Needle characteristics


- 4.8 Fr allows for use with smaller working channels and most pediatric scopes
- Black indicator tip provides visual proximity to tissue for accurate placement, safe penetration and gauged puncture depth
- Precise angle and 3-bevel tip engineered to minimize puncture force
- Unique cannula design provides easy injection and superior flexibility while maintaining structural integrity

## Available size options

- 35 cm – DIS199-E
- 70 cm – DIS201-E

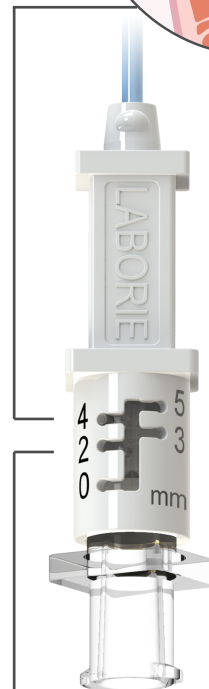

© 2024 Laborie. All rights reserved.

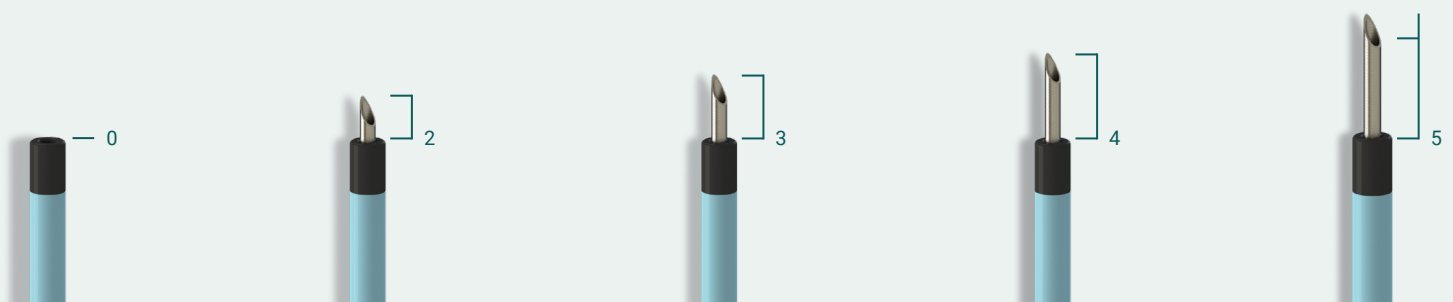

**Laborie**<sup>1</sup>

EUROPE/INTERNATIONAL

T +31 53 4803700  
E [info@laborie.com](mailto:info@laborie.com)  
W [laborie.com](http://laborie.com)

HEADQUARTERS

T +1 800 522 6742  
E [marketing@laborie.com](mailto:marketing@laborie.com)



Learn more